



Rest Assured with NiteOwl
Your companion for anytime
anywhere sleep therapy:

Camping

Maintain your sleep therapy
when camping or enjoying
other travel adventures.

*New swappable cell pack design
for limitless power options.*



*Powers your CPAP machine for 1-3 nights
with the heated humidifier turned off. *Run times depend on
many factors including device type, pressure setting, mask fit, etc.*

Air Travel

FAA-approved for use on board
an aircraft making it great for
international & red-eye flights.

*Small and lightweight so it's easy to
keep with your CPAP machine and supplies.*



3.0 USB PORT

Quick-charge 3.0 USB
port charges your go-to
electronic devices.

*Charge the NiteOwl with your
CPAP power supply - no need for extra cords.*



Dimensions:

8.7" L x 4.8" W x 0.9" H

Weight: 2 lbs.

Compatibility

NiteOwl is compatible with ResMed,
Philips Respironics, and ReactHealth's
Luna G3 and TravelPAP devices.



REACTHEALTH

NITEOWL PRODUCT CATALOG



866-718-0557



sleepniteowl@gmail.com



niteowlcpapbattery.com



NiteOwl Battery Kit

NO2-KIT

NiteOwl System Kit is compatible with the most popular CPAP units on the market. This Kit includes the (1) System Interface Assembly w/ cables, and (1) 150Wh Cell Pack:



ResMed Air 11



ResMed Air 10



Philips DreamStation 1



ResMed Mini



Philips DreamStation 2



Luna G3



Luna TravelPAP




NO2-150

NiteOwl 150Wh Cell Pack
Assembly Only

NO2-SIA

NiteOwl System Interface
Assembly Only



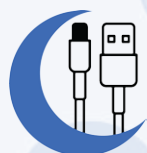

NiteOwl Connection Cables



PHRESP-CAB

NiteOwl connection cable for use with 12V DC Philips Respironics System One 60 series, DreamStation 1, and DreamStation 2 devices.

RESMED-CABKIT

Connection cables & power supply adapters for use with 24V DC ResMed Air 10, Air 11, and AirMini devices.

Kit includes the following cables:

- A10-CAB (Air 10 output cable)
- AM-A11-CAB (Air 11 and AirMini output cable)
- A11-PSA (Air 11 power supply adapter cable)
- AM-PSA (AirMini power supply adapter cable)



A10-CAB



AM-A11-CAB



A11-PSA



AM-PSA

LUNA-CABKIT

Luna G3 and Luna TravelPAP





NITE Owl
CPAP BATTERY

Rest Assured with NiteOwl
Your companion for anytime
anywhere sleep therapy:



NiteOwl Charging Accessories

AUTO-CAB

Mobile Power Adapter (car charger) to charge the NiteOwl via 12V DC outlet.



AC-CHARGER

100-240V AC wall charger when a separate power supply is preferred for charging the NiteOwl. *The NiteOwl is designed to use the PAP device AC power supply for charging.



NiteOwl Specifications

Swappable Cell Pack Design

Cell Type: Rechargeable Lithium-ion

Casing: Aluminum Alloy

Power Gauge: Digital Display (shows charge level as a %)

Capacity: 150 Wh

Output Voltage: 12V-24 DC (depending on DC output cable)

DC Input Voltage Range: 10V-28V

QC 3.0 USB Output Voltage: 5V | 2.4A (maximum)

Charge Input: 10V-28V

Optimal Operating Conditions: 23°F to 104°F

*capacity is reduced by 30% when operating between 23°F to 32°F

Optimal Charging Conditions: 32°F to 104°F

Storage Conditions: -4° to 104° (optimal range 60°F-80°F)

Charge Time: 2-3 hours @ 90W

*charge times are per battery using a 90W rated AC power supply

Life Span: 500 cycles (100%-0% then back to 100%)

*after 500 cycles the battery shall recover at least 70% of initial charge

Warranty: 1 year limited



866-718-0557



sleepniteowl@gmail.com



niteowlcpapbattery.com